



Forceps

Moorfields, non toothed/tying, notched, toothed, Capsulorhexis, Colibri, Cilia, Muscle/Strabismus, Lens, Artery, towel clamps, specialist



Forceps form an extension of the surgeons fingers enabling secure handling of tissue, suture, intraocular lens etc. They can be used to grasp, hold, pull/push, lift, or tear depending on the tip design and application.


Forceps anatomy

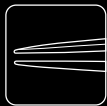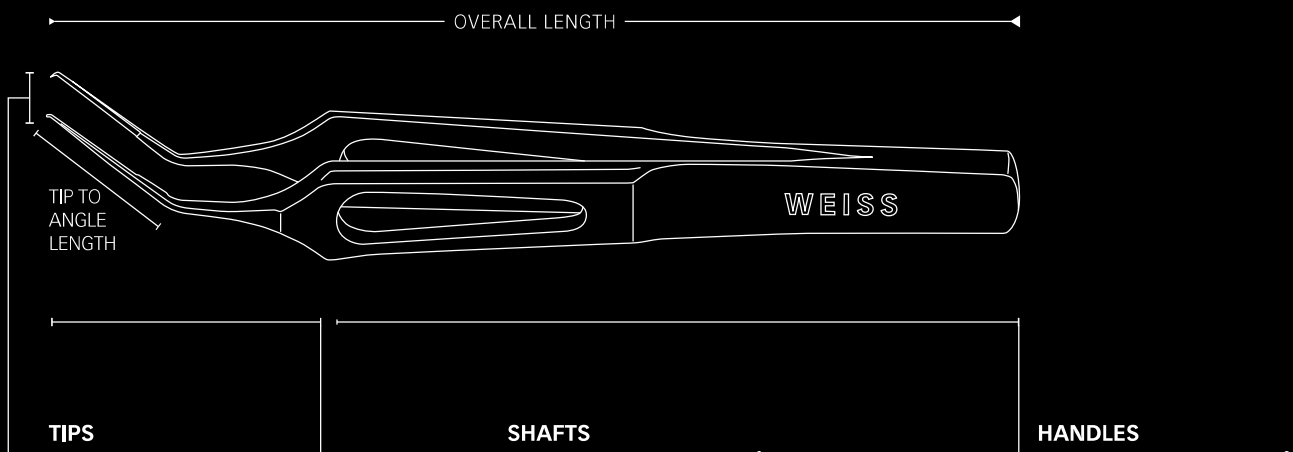


Plain –
general purpose –
removing sutures




Serrated –
general purpose –
grasping tissue (atraumatic)

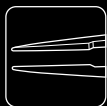


Platform –
for holding, tying and
removing sutures

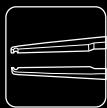

Notched –
for grasping and holding
tissue (atraumatic),
tying sutures (platform)

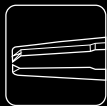

Toothed –
for grasping and holding
tissue including muscle,
tying sutures (platform)

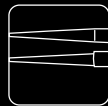

Straight –
for when access is direct




Curved –
for when access
is restricted




Angled –
for when access
is restricted




Colibri –
for good visualisation
and positioning of tip

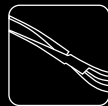

Cross action –
for maximum opening of tip
within small incision

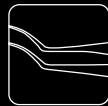

Formed –
accommodates natural
curvature of crystalline lens

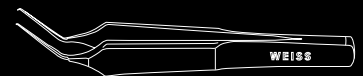

Flat – for good manual grip
and control of the tip but may
require significant hand, wrist
or forearm movement




Round – rotating handle by
rolling between fingers allows
for 'micro' control of the tip
with minimum movement



Cross action – opens and
closes tips by asserting and
releasing pressure – offering
'micro' control

Forceps – non toothed

Moorfields

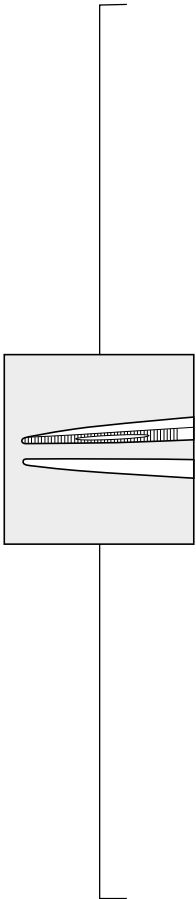

0101542 **SS**

Moorfields forceps, straight with standard serrated tips. 14mm platform, overall length 113mm



0109015 **Ti**

Moorfields forceps, straight with standard serrated tips. 14mm platform, overall length 113mm



0101543 **SS**

Moorfields forceps, straight with fine serrated tips. 14mm platform, overall length 113mm



0109021 **Ti**

Moorfields forceps, straight with fine serrated tips. 14mm platform, overall length 113mm



0101544 **SS**

Moorfields Skov forceps angled with standard tips. 14mm platform, 15mm tips to angle, overall length 110mm

Forceps – non toothed, tying

Plain

0109030 **Ti**

Weiss slimline forceps, straight with plain tips, overall length 85mm



0109024 **Ti**

Weiss forceps, straight with plain tips, overall length 85mm



0101166 **SS**

Hoskins forceps No.17, straight with plain 0.5mm tips, overall length 93mm



0101179 **SS**

Hoskins forceps No.32, straight with plain 0.3mm tips, overall length 88mm






0101181 **SS**

Hoskins No 34 Colibri forceps with fine pointed plain tips, no platform, overall length 86mm




0109428 **Ti**

Ring tipped forceps, atraumatic. Straight with 1.2mm ring and 0.6mm hole, overall length 108mm

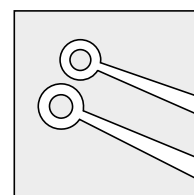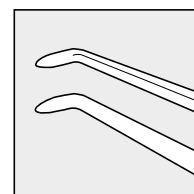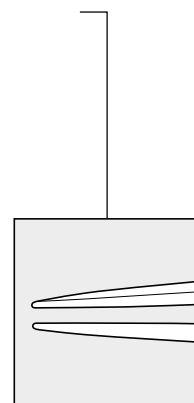

Forceps – non toothed, tying

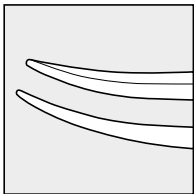


0101008 **SS**

Max Fine forceps, curved with plain tips, overall length 110mm



0109422 **Ti**

Max Fine forceps, curved with plain tips, overall length 110mm

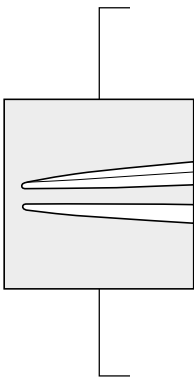

0101373 **SS**

Jewellers forceps No. 3, straight with plain, medium size tips, overall length 126mm




0101375 **SS**

Jewellers forceps No. 4, straight with plain, fine tips, overall length 106mm



0101372 **SS**

Jewellers forceps No. 5, straight with plain, ultra fine tips, overall length 115mm

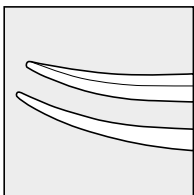

0101374 **SS**

Jewellers forceps No. 7, curved with plain tips, overall length 115mm

Discontinued

Discontinued

Discontinued

Forceps – non toothed, tying

Serrated

0101315 **SS**

Iris forceps, straight with 12mm serrated tips, overall length 99mm



0101323 **SS**

Graefe forceps straight with 12mm serrated tips, overall length 100mm



0204132 **SS**

Adson forceps, straight with 10mm serrated tips, overall length 122mm



0101324 **SS**

Graefe forceps curved, with 10mm serrated tips, overall length 93mm





0101313 **SS**

Iris forceps, 25° curved with 10mm serrated tips, overall length 99mm




0101314 **SS**

Iris forceps, 60° curved with 10mm serrated tips, overall length 99mm

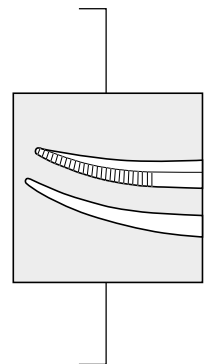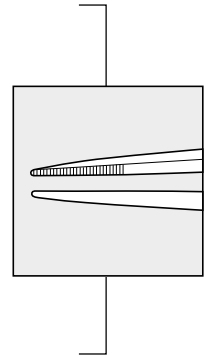

Forceps – non toothed, tying

Platform

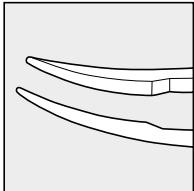

0109075 Ti


Weiss slimline forceps, curved with 4mm platform, overall length 83mm



0109073 Ti

Weiss slimline forceps, delicate tips, curved with 4mm platform, overall length 83mm

5mm Platform

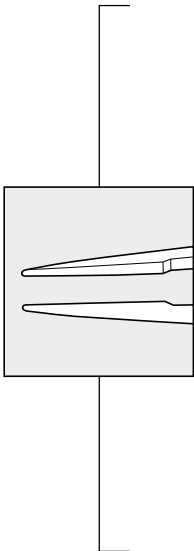

0101103 SS

Birks forceps, straight with 5mm platform, overall length 105mm



0109142 Ti

Birks forceps, straight with 5mm platform, overall length 109mm



0109035 Ti

Weiss slimline forceps, straight with 5mm platform, overall length 85mm



0109026 Ti

Weiss forceps, straight with 5mm platform, overall length 85mm


Forceps – non toothed, tying

Platform

5mm Platform continued

0101018 **SS**

Harms forceps, straight with 5mm platform, 0.5mm tips, overall length 101mm

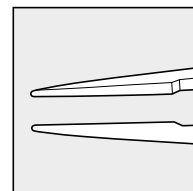


0101020 **SS**

Harms-Jaffe forceps, straight with 5mm platform, 0.3mm tips, overall length 101mm



0109064 **Ti**

Weiss Colibri forceps with 5mm platform, overall length 86mm

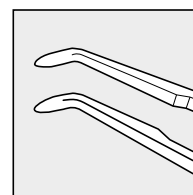


0101365 **SS**

Colibri forceps with 5mm platform, overall length 78mm



0101105 **SS**

Birks forceps, curved with 5mm platform, overall length 105mm

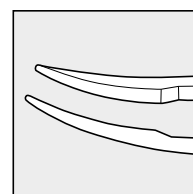

0109143 **Ti**



Birks forceps, curved with 5mm platform, overall length 109mm



Forceps – non toothed, tying

Platform

6mm Platform

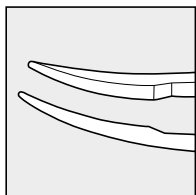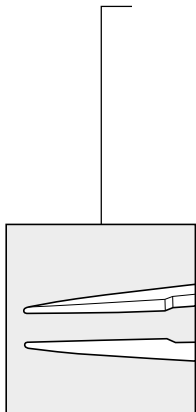

0101006 SS

Stallard forceps, straight with 6mm platform, overall length 111mm



0109152 Ti

Birks forceps, straight with 6mm platform, ultra fine tips, overall length 109mm



0101015 SS

Castroviejo wide handle forceps, straight with 6mm platform, delicate tips, overall length 105mm



0101010 SS

Catford forceps, straight with 6mm platform, overall length 73mm



0109420 Ti

Castroviejo forceps, straight with 6mm platform, delicate tips, overall length 110mm



0109421 Ti

Castroviejo forceps, curved with 6mm platform, delicate tips, overall length 110mm



0101019 SS

Harms forceps, curved with 6mm platform, overall length 100mm



0101021 SS

Harms-Jaffe forceps, curved with 6mm platform, overall length 100mm

Forceps – non toothed, tying

Platform

7mm Platform continued

0101017 **SS**

McPherson forceps, straight with 7mm platform, overall length 89mm



0101009 **SS**

Rycroft forceps, straight with 7mm platform, overall length 103mm



0101011 **SS**

McPherson forceps, angled with 7mm platform, 7.5mm tip to angle, overall length 91mm



0109097 **Ti**

McPherson forceps, angled with 7mm platform, 7.0mm tip to angle, overall length 85mm



0109149 **Ti**

Birks forceps, angled with 7mm platform, 7.5mm tip to angle, overall length 106mm



0109096 **Ti**

McPherson forceps, angled with 7mm platform, 10mm tip to angle, overall length 80mm



0109098 **Ti**

McPherson forceps, angled with 7mm platform, 12mm tip to angle, overall length 81mm



0109091 **Ti**



McPherson forceps, angled on flat with 7mm platform, 12mm tip to angle, overall length 85mm



NOTE: this forcep opens top to bottom



0109435 **Ti**

Kelman-McPherson forceps, angled with 7.5mm platform, 9mm tip to angle, overall length 108mm

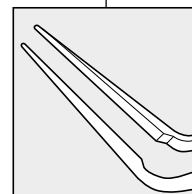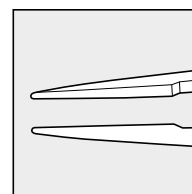

Forceps – non toothed, tying

Platform

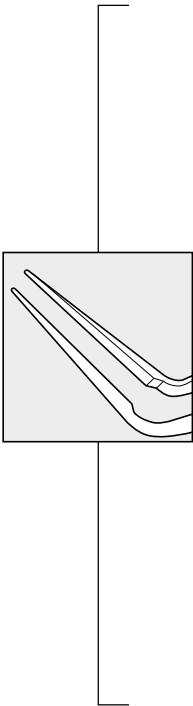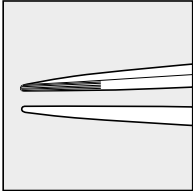

8mm Platform



0101003 SS

Alabama forceps, straight with 8mm platform and longitudinal serrations, overall length 101mm

10mm Platform



0101013 SS

Kelman-McPherson forceps, angled with 10mm platform, 10.5mm tip to angle, overall length 88mm



0109433 Ti

Kelman-McPherson forceps, angled with 10mm platform, 10.5mm tip to angle, overall length 103mm



0109148 Ti

Birks forceps, angled with 10mm platform, 10.5mm tip to angle, overall length 106mm

12mm Platform



0101014 SS

Kelman-McPherson forceps, angled with 12mm platform, 12.5mm tip to angle, overall length 88mm



0109434 Ti

Kelman-McPherson forceps, angled with 12mm platform, 12.5mm tip to angle, overall length 103mm

Forceps – notched

0109028 **Ti**

Weiss slimline forceps, straight with a 0.25mm notch, 4mm platform, overall length 85mm



0109023 **Ti**

Weiss forceps, straight with a 0.25mm notch, 3mm platform, overall length 85mm



0109427 **Ti**

Castroviejo forceps, straight with a 0.25mm notch, 6mm platform, overall length 110mm



0109072 **Ti**

Weiss forceps, curved with a 0.25mm notch, 4mm platform, overall length 83mm



0109150 **Ti**

Birks forceps, angled with a 0.25mm notch, 7mm platform, 8mm tip to angle, overall length 106mm






0101012 **SS**

McPherson forceps, angled with a 0.3mm notch, 7mm platform, 7.5mm tip to angle, overall length 91mm





0109431 **Ti**

McPherson forceps, angled with a 0.3mm notch, 7mm platform, 7mm tip to angle, overall length 108mm

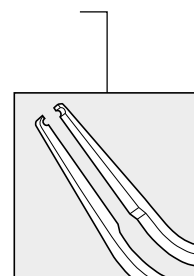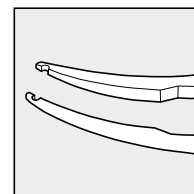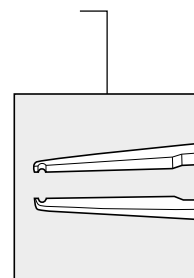

Forceps – notched

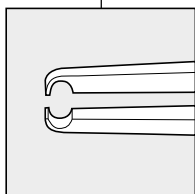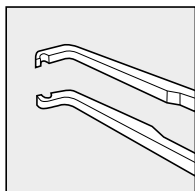

0109147 **Ti**

Birks Colibri forceps, with a 0.25mm notch, 4mm platform, overall length 108mm



0109063 **Ti**

Weiss Colibri Forceps, with a 0.25mm notch, 5mm platform, overall length 86mm



0101162 **SS**

Hoskins No.10 forceps, straight with a 0.8mm notch, overall length 93mm



0101163 **SS**

Hoskins No.11 forceps, straight with a 0.5mm notch, overall length 93mm



0101164 **SS**

Hoskins No. 15 forceps, straight with a 0.3mm notch, overall length 93mm



0101176 **SS**

Hoskins No.28 forceps, straight with a 0.25mm notch, overall length 88mm



0101168 **SS**

Hoskins No.19 forceps, straight with a 0.3mm notch, overall length 88mm



0101167 **SS**

Hoskins No. 18 forceps, straight with a 0.5mm notch, overall length 88mm



Forceps – toothed

0.1mm teeth

0109036 Ti
Weiss slimline forceps, straight with 0.1mm forward facing teeth, 5mm platform, overall length 85mm



0109065 Ti
Weiss Colibri forceps, with 0.1mm forward facing teeth, 5mm platform, overall length 85mm

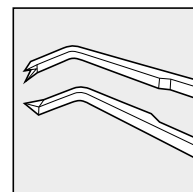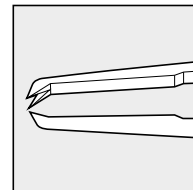

0.12mm teeth

0101107 SS
Birks forceps, straight with 0.12mm teeth, 5mm platform, overall length 105mm



0109144 Ti
Birks forceps, straight with 0.12mm teeth, 5mm platform, overall length 109mm



0101075 SS
Bonn forceps, straight with 0.12mm teeth, 6mm platform, overall length 115mm



0109013 Ti
Bonn forceps, straight with 0.12mm teeth, 6mm platform, overall length 111mm





0109012 Ti
Jayle forceps, straight with 0.12mm teeth, 6mm platform, overall length 118mm



0109014 Ti
Bonn forceps, straight with 0.12mm teeth, overall length 72mm






0101332 SS
Bonn forceps straight with 0.12mm teeth, overall length 72mm

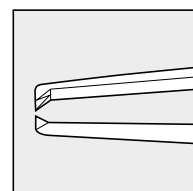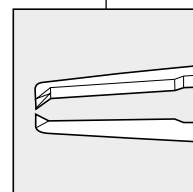

Forceps – toothed

0.12mm teeth

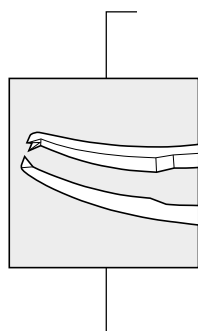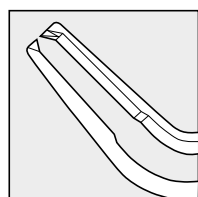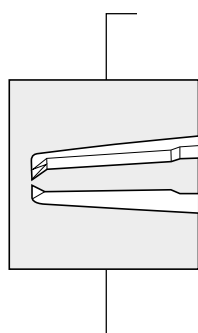

0109029 Ti

Weiss slimline forceps, straight with 0.12mm teeth, 4mm platform, overall length 85mm



0109025 Ti

Weiss forceps, straight with 0.12mm teeth, 5mm platform, overall length 85mm



0109424 Ti

Castroviejo forceps, straight with 0.12mm teeth, 6mm platform, overall length 110mm



0109095 Ti

McPherson forceps, angled with 0.12mm teeth, 7mm platform, 14mm tip to angle, overall length 83mm



0109151 Ti

Birks forceps, angled with 0.12mm teeth, 10mm platform, 10.5mm tip to angle, overall length 106mm



0101109 SS

Birks forceps, curved with 0.12mm teeth, 5mm platform, overall length 105mm



0109145 Ti

Birks forceps, curved with 0.12mm teeth, 5mm platform, overall length 109mm



0109074 Ti

Weiss forceps, curved with 0.12mm teeth, 5mm platform, overall length 83mm

Forceps – toothed

0.12mm teeth continued

0101111 SS

Birks Colibri forceps with 0.12mm teeth, 5mm platform, overall length 105mm



0109146 Ti

Birks Colibri forceps with 0.12mm teeth, 4mm platform, overall length 107mm



0101079 SS

Castroviejo Colibri forceps with 0.12mm teeth, 5mm platform, overall length 105mm



0109423 Ti

Castroviejo Colibri wide handle forceps with 0.12mm teeth, 5mm platform, overall length 110mm



0101364 SS

Colibri forceps with 0.12mm teeth, 6mm platform, overall length 74mm



0109062 Ti

Weiss Colibri forceps with 0.12mm teeth, 5mm platform, overall length 86mm



0101074 SS

Castroviejo forceps, straight with 0.12mm forward facing teeth, 8mm platform, overall length 113mm






0101078 SS

Castroviejo wide handle forceps, straight with 0.12mm forward facing teeth, 6mm platform, overall length 105mm






0109070 Ti

Polack Colibri forceps, 2 pronged for corneal suturing with 0.12mm teeth, 4.5mm platform, overall length 75mm

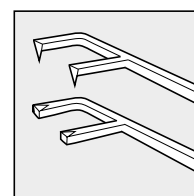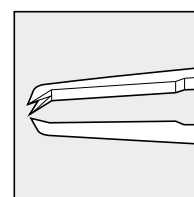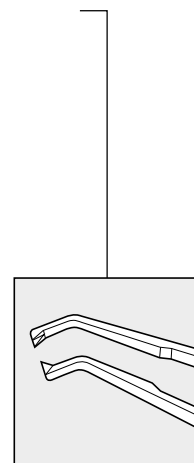

Forceps – toothed

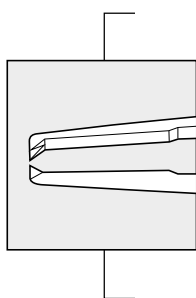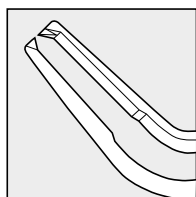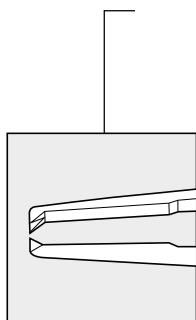

0101072 SS

Jayle forceps, straight with 0.15mm teeth, 6mm platform, overall length 118mm



0101073 SS

St. Martin forceps, straight with 0.15mm teeth, 6mm platform, overall length 87mm



0109162 Ti

St. Martin forceps, straight with 0.15mm teeth, 5mm platform, overall length 88mm



0101076 SS

Catford forceps, straight with 0.15mm teeth, 4mm platform, overall length 70mm



0101077 SS

McPherson forceps, angled with 0.15mm teeth, 7mm platform, 7.5mm tip to angle, overall length 91mm

0.3mm teeth



0101045 SS

Rycroft forceps, straight with 0.3mm teeth, 8mm platform, overall length 104mm



0101042 SS

Jayle forceps, straight with 0.3mm teeth, 6mm platform, overall length 118mm



0109011 Ti

Jayle forceps, straight with 0.3mm teeth, 6mm platform, overall length 115mm

Forceps – toothed

0.3mm teeth continued

0101043 SS

St.Martin forceps, straight with 0.3mm teeth, 6mm platform, overall length 87mm



0109172 Ti

St.Martin forceps, straight with 0.3mm teeth, 5mm platform, overall length 88mm



0101049 SS

Castroviejo forceps, straight with 0.3mm teeth, 6mm platform, overall length 99mm



0109425 Ti

Castroviejo forceps, straight with 0.3mm teeth, 6mm platform, overall length 110mm



0101050 SS

Castroviejo forceps, straight with 0.3mm teeth, 6mm platform, overall length 99mm



0109180 Ti

Birks forceps, curved with 0.3mm teeth, 5mm platform, overall length 109mm



0101362 SS

Colibri forceps with 0.3mm teeth, overall length 74mm







0101363 SS

Colibri forceps with 0.3mm teeth, 5mm platform, overall length 74mm




0101046 SS

Castroviejo wide handle forceps, straight with 0.3mm forward facing teeth, 8mm platform, overall length 105mm

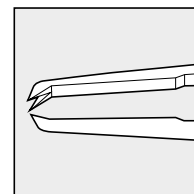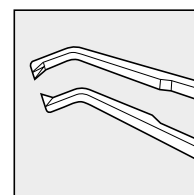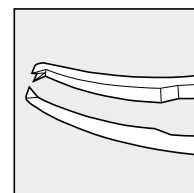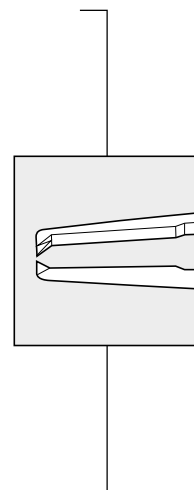


Forceps – toothed

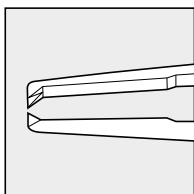

0.5mm teeth


0109426 **Ti**

Castroviejo forceps, straight with 0.5mm teeth, 6mm platform, overall length 110mm

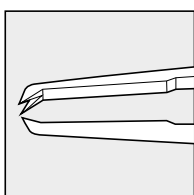


0101047 **SS**

Castroviejo wide handle forceps, straight with 0.5mm forward facing teeth, 8mm platform, overall length 105mm

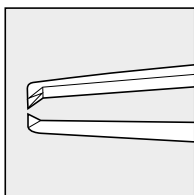

0204134 **SS**


Adson forceps, straight with 0.5mm teeth, overall length 122mm

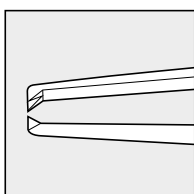

0.75mm teeth

0101318 **SS**

Graefe forceps, straight with 0.75mm teeth, overall length 73mm

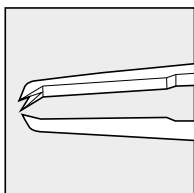

0.9mm teeth

0101048 **SS**


Castroviejo wide handle forceps, straight with 0.9mm forward facing teeth, 8mm platform, overall length 105mm

Forceps – toothed

Fixation

0101284 **SS**

Weiss fixation forceps, straight with 1mm teeth (1 into 2), overall length 115mm

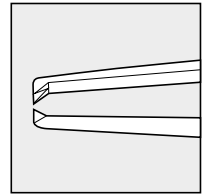

0101285 **SS**

Francis fixation forceps, straight with 1mm teeth (2 into 3), overall length 120mm




0101298 **SS**

Graefe locking fixation forceps, straight with 4mm wide tip, 6 x 0.7mm teeth, overall length 110mm




0101005 **SS**

Lister forceps, straight with 0.65mm forward facing teeth, 12mm platform, overall length 110mm

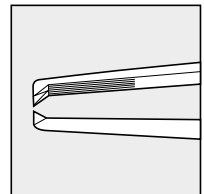


0101290 **SS**

Elschnig forceps, straight with 1.2mm teeth (1 into 2), overall length 109mm

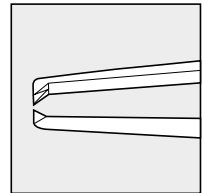


Forceps – Capsulorhexis

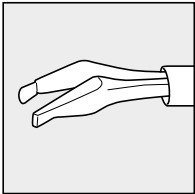

For incision size 1mm - 1.8mm


2003100 Ti

23G, curved, capsulorhexis forceps, end grasping for sub 1mm incision, overall length 140mm

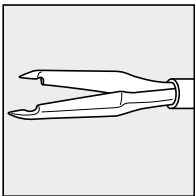

2003102 Ti

Gutierrez-C combination curved 23G capsulorhexis forceps/scissors for sub 1mm incision. 0.3mm end grasping forcep lies in front of straight scissor with 2.0mm vertical cutting blade, overall length 140mm

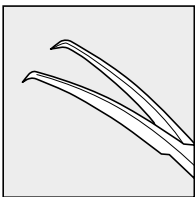

For incision size 1.8mm - 2.0mm


0109081 Ti

Cross action capsulorhexis forceps with sharp tips. Formed shaft, 14mm tip to angle, overall length 105mm



0109083 Ti

Cross action round body forceps with sharp tips. Formed shaft, 14mm tip to angle, overall length 115mm

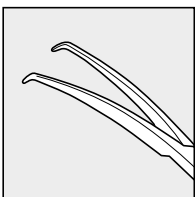

0109082 Ti

Cross action capsulorhexis forceps with round tips. Formed shaft, 14mm tip to angle, overall length 105mm

Forceps – Capsulorhexis

For incision size 2mm - 2.5mm

0101569 **SS**

Capsulorhexis forceps, round body with ultrafine tips. Formed shaft, 9mm tip to angle, overall length 104mm




0101570 **SS**

Capsulorhexis forceps, wide body with ultrafine tips. Formed shaft, 10mm tip to angle, overall length 105mm



0109136 **Ti**

Weiss capsulorhexis forceps with ultrafine, serrated tips and 1mm markings along shaft. Formed shaft, 14mm tip to angle, overall length 80mm

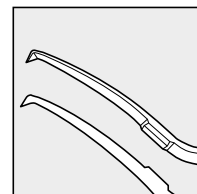

For incision size 2.5mm - 2.8mm


0109127 **Ti**

Inamura round bodied, cross action capsulorhexis forceps with fine tapered tips. Formed shaft, 14mm tip to angle, overall length 110mm



0109129 **Ti**

Inamura round bodied, cross action capsulorhexis forceps with fine tapered tips. Formed shaft, 13mm tip to angle, overall length 110mm

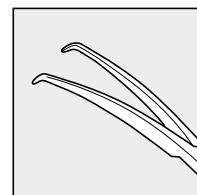

For incision size 2.8mm+ - with formed shaft


0101563 **SS**

Capsulorhexis forceps, round body with standard tips. Formed shaft, 12mm tip to angle, overall length 103mm




0109122 **Ti**

Capsulorhexis forceps, round body with fine tips. Formed shaft, 12mm tip to angle, overall length 105mm

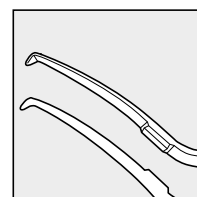

Forceps – Capsulorhexis

For incision size 2.8mm+ - with formed shaft

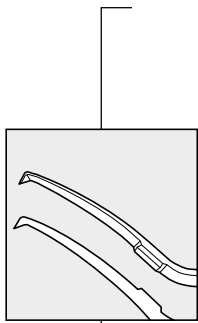

0109123 **Ti**

Capsulorhexis forceps, round body with sharp tips. Formed shaft, 12mm tip to angle, overall length 105mm



0109445 **Ti**

Utrata capsulorhexis forceps with sharp tips. Formed shaft, 12mm tip to angle, overall length 104mm



0109117 **Ti**

Weiss capsulorhexis forceps with sharp tips. Formed shaft, 12mm tip to angle, overall length 108mm



0109119 **Ti**

Cross action capsulorhexis forceps with sharp tips. Formed shaft, 12mm tip to angle, overall length 108mm



0109438 **Ti**

Utrata capsulorhexis forceps with fine tips. Formed shaft, 12mm tip to angle, overall length 104mm



0101565 **SS**

Capsulorhexis forceps, wide body with standard tips. Formed shaft, 12mm tip to angle, overall length 105mm




0109113 **Ti**

Weiss capsulorhexis forceps with fine tips. Formed shaft, 12mm tip to angle, overall length 108mm



0109115 **Ti**

Cross action capsulorhexis forceps with fine tips. Formed shaft, 12mm tip to angle, overall length 108mm

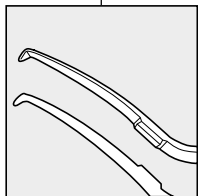

Forceps – Capsulorhexis

For incision size 2.8mm+ - with standard shaft

0109135 **Ti**

Weiss slimline capsulorhexis forceps with sharp, serrated tips. Standard shaft, 14mm tip to angle, overall length 80mm



0109444 **Ti**

Utrata capsulorhexis forceps with sharp tips. Standard shaft, 12mm tip to angle, overall length 104mm



0109120 **Ti**

Capsulorhexis forceps, round body with sharp tips. Standard shaft, 12mm tip to angle, overall length 105mm



0101562 **SS**

Capsulorhexis forceps, round body with standard tips. Standard shaft, 12mm tip to angle, overall length 103mm



0109121 **Ti**

Capsulorhexis forceps, round body with fine tips. Standard shaft, 12mm tip to angle, overall length 103mm





0109131 **Ti**

Weiss slimline capsulorhexis forceps with standard, serrated tips. Standard shaft, 14mm tip to angle, overall length 80mm




0109437 **Ti**

Utrata capsulorhexis forceps with fine tips. Standard shaft, 12mm tip to angle, overall length 104mm

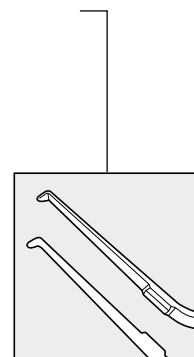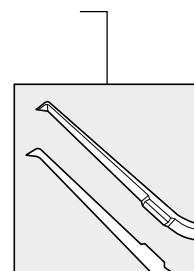

Forceps – Capsulorhexis

For incision size 2.8mm+ - with standard shaft

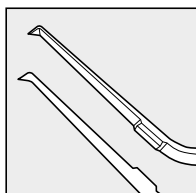

0109116 Ti

Weiss capsulorhexis forceps with sharp tips. Standard shaft, 12mm tip to angle, overall length 108mm



0109118 Ti

Cross action capsulorhexis forceps with sharp tips. Standard shaft, 12mm tip to angle, overall length 108mm



0101564 SS

Capsulorhexis forceps, wide body with standard tips. Standard shaft, 12mm tip to angle, overall length 106mm




0109112 Ti

Weiss capsulorhexis forceps with fine tips. Standard shaft, 12mm tip to angle, overall length 108mm



0109114 Ti


Cross action capsulorhexis forceps with fine tips. Standard shaft, 12mm tip to angle, overall length 108mm

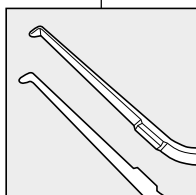

Forceps – Colibri


0101181 **SS**

Hoskins No 34 Colibri forceps with fine pointed plain tips, no platform, overall length 86mm

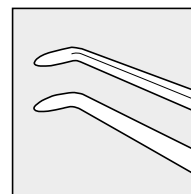


0101362 **SS**

Colibri forceps with 0.3mm teeth, no platform, overall length 74mm

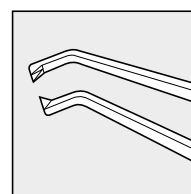


0109147 **Ti**

Birks Colibri forceps, with a 0.25mm notch, 4mm platform, overall length 108mm

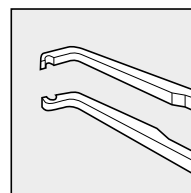


0109063 **Ti**

Weiss Colibri forceps, with a 0.25mm notch, 5mm platform, overall length 86mm

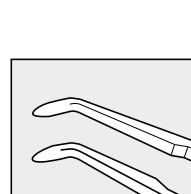


0109064 **Ti**

Weiss Colibri forceps with 5mm platform, overall length 86mm

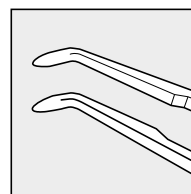


0101365 **SS**

Colibri forceps with 5mm platform, overall length 78mm

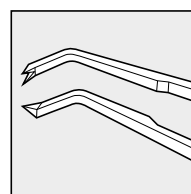

0109065 **Ti**

Weiss Colibri forceps, with 0.1mm forward facing teeth, 5mm platform, overall length 85mm

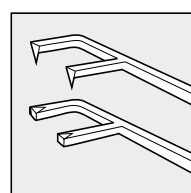


0109070 **Ti**

Polack Colibri forceps, 2 pronged for corneal suturing with 0.12mm teeth, 4.5mm platform, overall length 75mm



Forceps – Colibri

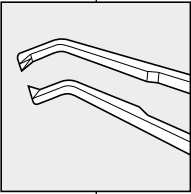

0101079 SS

Castroviejo Colibri wide handle forceps with 0.12mm teeth, 5mm platform, overall length 105mm



0109423 Ti

Castroviejo Colibri wide handle forceps with 0.12mm teeth, 5mm platform, overall length 110mm



0101364 SS

Colibri forceps with 0.12mm teeth, 6mm platform, overall length 74mm



0109062 Ti

Weiss Colibri forceps with 0.12mm teeth, 5mm platform, overall length 86mm



0101363 SS

Colibri forceps with 0.3mm teeth, 5mm platform, overall length 74mm



0101111 SS

Birks Colibri forceps with 0.12mm teeth, 5mm platform, overall length 105mm



0109146 Ti

Birks Colibri forceps with 0.12mm teeth, 4mm platform, overall length 107mm

Cilia Forceps

0101264 **SS**

Whitfield cilia forceps, straight with angled tip, 3mm platform, overall length 91mm



0101265 **SS**

Bennett cilia forceps, straight with 3mm cup tips, overall length 75mm



0101266 **SS**

Cilia forceps, straight with rounded tips, 4mm platform, overall length 78mm



0101268 **SS**

Cross action cilia forceps, angled with rounded tips, 11mm tip to angle, overall length 105mm



0101267 **SS**

Barraquer cilia forceps, angled with rounded tips, 7mm platform, 7mm tip to angle, overall length 85mm



0101269 **SS**

Barraquer cilia forceps, angled with with rounded tips, 8mm platform, 8mm tip to angle, overall length 100mm



Muscle/Strabismus Forceps



0101592 SS

Troutman superior rectus forceps, angled with 0.5mm teeth. 10mm tip to angle, overall length 105mm



0101291 SS

Dastoor rectus forceps, 'S' shaped with 1 into 2 teeth. 10mm tip to angle, overall length 110mm



0101396 SS

Jameson strabismus forceps, locking, adult, angled **right** with 6 x 0.6mm teeth. 13mm tip to angle, overall length 99mm

0101394 SS

Jameson strabismus forceps, locking, adult, angled **left** with 6 x 0.6mm teeth. 13mm tip to angle, overall length 99mm



0101397 SS

Jameson strabismus forceps, locking, child, angled **right** with 4 x 0.6mm teeth. 8mm tip to angle, overall length 99mm

0101395 SS

Jameson strabismus forceps, locking, child, angled **left** with 4 x 0.6mm teeth. 18mm tip to angle, overall length 99mm



0101398 SS

Worth strabismus forceps, locking, angled **right** with 15mm longitudinal serrations. 13mm tip to angle, overall length 100mm

0101399 SS

Worth strabismus forceps, locking, angled **left** with 15mm longitudinal serrations. 13mm tip to angle, overall length 100mm




0109531 Ti

Federio müllerectomy forceps, 'T' ended with 25mm bar, 12 x 0.8mm interlocking teeth, overall length 105mm

Lens Forceps


0101593 **SS**
 IOL folding forceps, angled with grooved paddles.
 10mm tip to angle, overall length 108mm

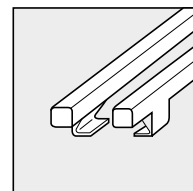

0109440 **Ti**
 IOL folding forceps, angled with grooved paddles
 and central marker. 8mm tip to angle, overall length
 108mm




0109163 **Ti**
 IOL holding forceps, cross action, notch tip, angled
 with a catch. 25mm tip to angle, overall length
 118mm

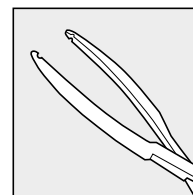


0109164 **Ti**
 IOL holding forceps, cross action, notch tip, angled,
 no catch. 25mm tip to angle, overall length 118mm



0109448 **Ti**
 IOL loading forceps, angled with atraumatic
 rounded tips and jaws, 15mm tip to angle, overall
 length 105mm

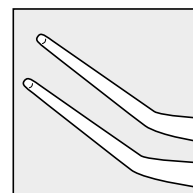

0109449 **Ti**
 IOL loading forceps, straight with atraumatic rounded
 tips and jaws, 15mm jaws, overall length 110mm

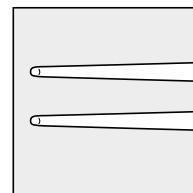

Artery Forceps



0101472 SS

Halstead mosquito forceps, straight with 20mm serrated tips, overall length 125mm



0101473 SS

Halstead mosquito forceps, curved with 26mm serrated tips, overall length 125mm



0101474 SS

Wells artery forceps, straight with 20mm serrated tips, overall length 95mm



0101475 SS

Wells artery forceps, curved with 20mm serrated tips, overall length 95mm

Forceps – towel clamps

0101589 **SS**

Towel clamp, cross action, straight, overall length 85mm



0101583 **SS**

Towel clamp, non-penetrating, curved, 5mm ball & socket tip, overall length 90mm



0204136 **SS**

Ramsey's sponge forceps, straight, 26mm serrated tips, overall length 182mm




0204147 **SS**

Backhaus towel clamp, curved, overall length 89mm



Forceps – specialist

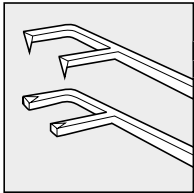

0109070 **Ti**

Polack Colibri forceps, 2 pronged for corneal suturing with 0.12mm teeth, 4.5mm platform, overall length 75mm



0101605 **SS**

Scleral plug forceps, cross action, angled tips have a 2mm groove for gripping scleral plug. 8mm tip to angle, overall length 106mm



0101552 **SS**

Watski sleeve forceps, reverse cross action for applying sleeve to retinal bands. Curved tips, overall length 116mm



0109181 **Ti**

Felipe fixation forceps, reverse cross action for fixating the eye during phaco. Angled, 12mm tip to angle, overall length 106mm



0101462 **SS**

Serrafine cross action mini clamp, straight. 10mm serrated jaws, overall length 38mm



0101463 **SS**

Serrafine cross action, feather weight, mini clamp, straight. 8mm serrated jaws, overall length 32mm



0101464 **SS**

Serrafine cross action, mini clamp, curved. 16mm serrated jaws, overall length 48mm

Specialist

0109185 **Ti**

Khaw conjunctival 'T' ended clamp with smooth, 12mm bar, overall length 86mm



0109186 **Ti**

Khaw-Shah conjunctival 'T' ended clamp with serrated 4mm bar, overall length 86mm



0204132 **SS**

Adson forceps, straight with 10mm serrated tips, overall length 122mm



0204134 **SS**

Adson forceps, straight with 0,5mm teeth, overall length 122mm

