




Scissors

Spring, bow/ ribbon. De Wecker, Vannas, Westcott, Williamson-Noble, Corneal, Capsulotomy, Strabismus, Stevens tenotomy, specialist, utility

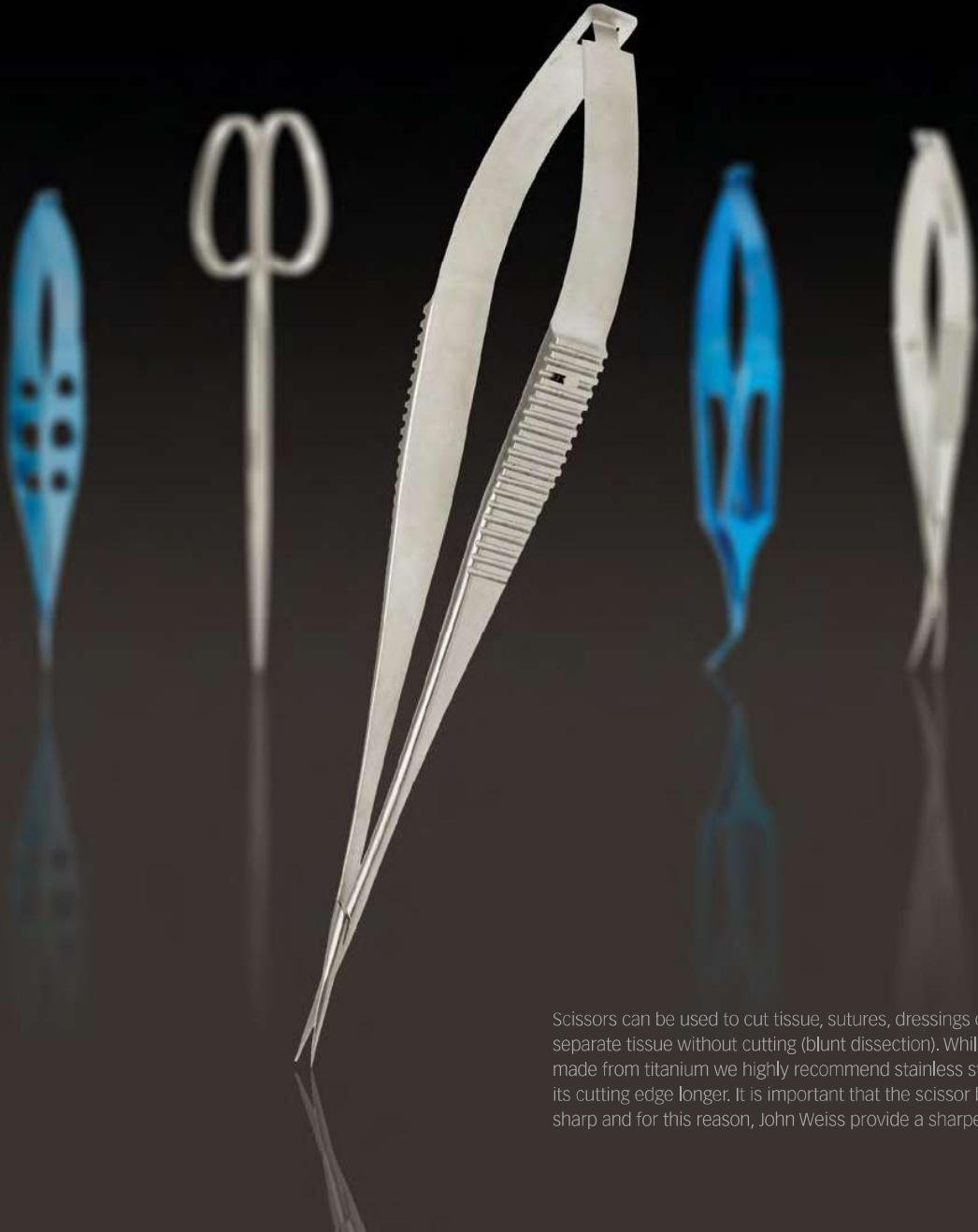


Scissors can be used to cut tissue, sutures, dressings or to divide or separate tissue without cutting (blunt dissection). While scissors can be made from titanium we highly recommend stainless steel as it maintains its cutting edge longer. It is important that the scissor blades are kept sharp and for this reason, John Weiss provide a sharpening service.

Scissors anatomy

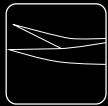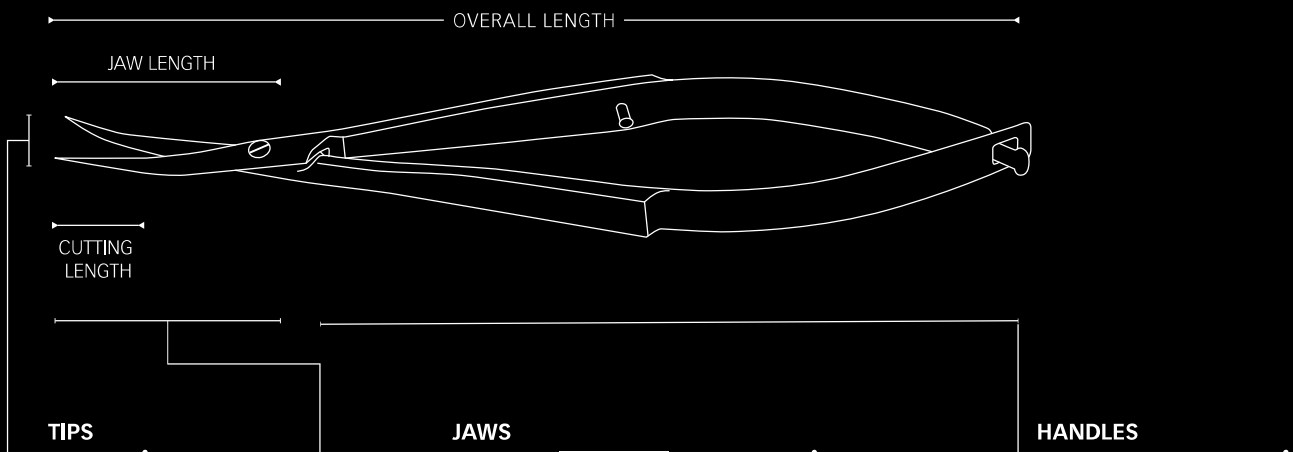

Sharp/pointed – can pierce tissue as well as cut

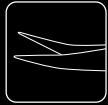


Rounded – for cutting only




Straight – good for when direct access is available or when a straight cut is required



Curved – to reach more difficult areas or when a curved cut is required. May allow better visualisation of tips

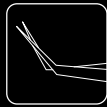

Angled – enables blades to be positioned at the correct angle to the tissue, particularly useful when excising

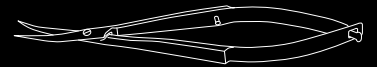

Spring – a good manual grip with excellent control of the tip means that these are the best choice for delicate work.



Bow or ribbon – traditional design, good utility scissor or for non-delicate excisions such as extraocular. May require significant hand and wrist movement.



Hinge – a time tested design, blades close by squeezing handle together. Generally used for iris excision.

Spring Scissors

De Wecker



0103245 SS

De Wecker Iris Scissors. Angled on flat with sharp tips. Max cutting length 5mm, blade length 8mm, overall length 110mm



0103242 SS

De Wecker iris scissors. Angled on flat with sharp tips. Max cutting length 9mm, blade length 12mm, overall length 115mm



0103243 SS

De Wecker Iris Scissors. Angled on flat with right tip sharp and left tip rounded. Max cutting length 9mm, blade length 12mm, overall length 115mm



0103244 SS

De Wecker Iris Scissors. Angled on flat with left tip sharp and right tip rounded. Max cutting length 9mm, blade length 12mm, overall length 115mm

0103252 SS

Barraquer Iris Scissors. Angled on flat with sharp tips. Max cutting length 6mm, blade length 8.5mm, overall length 58mm

Spring Scissors

Vannas – straight

0103102 **SS**

Vannas Scissors, straight with sharp tips.
Max cutting length 4mm, blade length 9mm,
overall length 86mm



0103105 **SS**

Vannas Scissors, straight with sharp tips.
Max cutting length 6mm, blade length 12mm,
overall length 89mm



0109003 **Ti**

Vannas Scissors, straight with sharp tips.
Max cutting length 3mm, blade length 9mm,
overall length 95mm



0109460 **Ti**

Vannas Scissors, straight with ultrafine, sharp
tips. Max cutting length 3mm, blade length 8mm,
overall length 110mm



0103123 **SS**

Miniature Scissors, straight with sharp tips.
Max cutting length 3mm, blade length 6.5mm,
overall length 85mm



0103129 **SS**

Miniature Scissors, straight sharp tips, long
handled. Max cutting length 3mm, blade
length 9mm, overall length 110mm



0103143 **SS**

Miniature Scissors, straight sharp tips, long
handled. Max cutting length 5mm, blade length
12mm, overall length 115mm



Spring Scissors

Vannas – curved



0103103 **SS**

Vannas Scissors, curved with sharp tips.
Max cutting length 4mm, blade length 9mm,
overall length 87mm



0103106 **SS**

Vannas Scissors, curved with sharp tips.
Max cutting length 5mm, blade length 12mm,
overall length 89mm



0109005 **Ti**

Vannas Scissors, curved with sharp tips.
Max cutting length 4mm, blade length 9mm,
overall length 94mm



0109461 **Ti**

Vannas Scissors, curved with ultrafine, sharp
tips. Max cutting length 3mm, blade length 6mm,
overall length 106mm



0109462 **Ti**

Vannas Scissors, curved with ultrafine sharp tips
angled forward. Max cutting length 3mm, blade
length 6mm, overall length 106mm



0103128 **SS**

Miniature Scissors, curved with sharp tips.
Max cutting length 3mm, blade length 6.5mm,
overall length 82mm



0103122 **SS**

Miniature Scissors, angled with sharp tips,
angled on flat. Max cutting length 3mm, blade
length 9mm, overall length 82mm



0103124 **SS**

Miniature Scissors, angled with sharp tips side
cutting. Max cutting length 3mm blade length
9mm, overall length 82mm

Spring Scissors

Westcott

0103172 SS

Westcott Scissors, curved with rounded tips. Max cutting length 9mm, blade length 23mm, overall length 115mm



0103130 SS

McPherson-Westcott Scissors, curved with sharp tips. Max cutting length 6mm, blade length 12mm, overall length 100mm



0103173 SS

Westcott Scissors, curved with sharp tips. Max cutting length 9mm, blade length 23mm, overall length 115mm



0103182 SS

Westcott Scissors, curved with rounded tips. Wide handle with a max cutting length 12mm, blade length 23mm, overall length 139mm



0103183 SS

Westcott Scissors, curved with sharp tips. Wide handle with a max cutting length 12mm, blade length 23mm, overall length 139mm



Spring Scissors

Williamson-Noble



0103160 SS

Williamson-Noble Strabismus Scissors. Straight with rounded tips. Max cutting length 10mm, blade length 27mm, overall length 125mm



0103161 SS

Williamson-Noble Strabismus Scissors. Curved with rounded tips. Max cutting length 10mm, blade length 27mm, overall length 125mm



0103230 SS

Williamson-Noble Iris Scissors. Straight with sharp tips. Max cutting length 14mm, blade length 27mm, overall length 115mm



0103231 SS

Williamson-Noble Iris Scissors. Curved with sharp tips. Max cutting length 14mm, blade length 27mm, overall length 115mm

Spring Scissors

Corneal

0103033 **SS**

Castroviejo Corneal Scissors curved with fine rounded tips. Max cutting length 6mm, blade length 12mm, overall length 107mm



0103032 **SS**

Castroviejo Corneal Scissors curved with rounded tips. Max cutting length 9mm, blade length 21mm, overall length 120mm



0103040 **SS**

Castroviejo Corneal Scissors. Miniature, curved, wide handle with sharp tips. Max cutting length 5mm, blade length 14mm, overall length 129mm



Right or left cutting

0103019 **SS**

Katzin Corneal Scissors. Angled right cutting with fine rounded tips. Max cutting length 5mm, blade length 10mm, overall length 103mm



0103020 **SS**

Katzin Corneal Scissors. Angled left cutting with fine rounded tips. Max cutting length 5mm, blade length 10mm, overall length 103mm



0103046 **SS**

Troutman-Katzin Corneal Scissors Angled right cutting with rounded tips. Max cutting length 4mm, blade length 9mm, overall length 100mm



0103047 **SS**

Troutman-Katzin Corneal Scissors Angled left cutting with rounded tips. Max cutting length 4mm, blade length 9mm, overall length 100mm



Spring Scissors



0103012 **SS**

Castroviejo Corneal Scissors Angled right cutting with rounded tips. Max cutting length 8mm, blade length 15mm, overall length 107mm



0103015 **SS**

Castroviejo Corneal Scissors Angled left cutting with rounded tips. Max cutting length 8mm, blade length 15mm, overall length 107mm

Capsulotomy



0103082 **SS**

Capsulotomy scissors. Curved with sharp tips. Max cutting length 6mm, blade length 10.5mm, overall length 105mm



0103083 **SS**

Capsulotomy scissors round body. Curved with sharp tips. Max cutting length 6mm, blade length 12mm, overall length 122mm



0103085 **SS**

Ong scissors. Angled on flat with sharp tips. Max cutting length 7mm, blade length 14mm, overall length 90mm



0109002 **Ti**

Ong scissors. Angled on flat with sharp tips. Max cutting length 8mm, blade length 10mm, overall length 95mm



0103088 **SS**

Gills-Welsh scissors, Angled on flat with sharp tips Max cutting length 6mm, blade length 13mm, overall length 84mm

Bow / Ribbon Style Scissors

Strabismus

0103152 **SS**

Strabismus Scissors, straight with rounded tips. Max cutting length 25mm, blade length 32mm, overall length 115mm



0103162 **SS**

Strabismus Scissors, short shank, straight with rounded tips. Max cutting length 19mm, blade length 28mm, overall length 89mm



0103156 **SS**

Strabismus Scissors, ribbon style, straight with rounded tips. Max cutting length 24mm, blade length 31mm, overall length 110mm



0103159 **SS**

Strabismus Scissors, ribbon style, curved with ultrafine rounded tips. Max cutting length 17mm, blade length 19mm, overall length 98mm



0103154 **SS**

Strabismus Scissors, curved with rounded tips. Max cutting length 20mm, blade length 26mm, overall length 90mm



Bow / Ribbon Style Scissors

Stevens Tenotomy



0103174 **SS**

Stevens Tenotomy Scissors, straight with rounded tips. Max cutting length 16mm, blade length 22mm, overall length 104mm



0103175 **SS**

Stevens Tenotomy Scissors, curved with rounded tips. Max cutting length 19mm, blade length 22mm, overall length 105mm



0103179 **SS**

Stevens Tenotomy Scissors, curved with rounded tips. Max cutting length 24mm, blade length 31mm, overall length 112mm



0103178 **SS**

Stevens Tenotomy Scissors, straight with rounded tips. Max cutting length 26mm, blade length 31mm, overall length 115mm

Bow / Ribbon Style Scissors

Special

0103002 **SS**

Pooley Conjunctival Scissors, straight, rounded tips. Max cutting length 18mm, blade length 26mm, overall length 90mm



0103003 **SS**

Pooley Conjunctival Scissors, straight, rounded tips. Max cutting length 22mm, blade length 30mm, overall length 113mm



0103222 **SS**

Landolt Iris Scissors, curved with rounded tips. Max cutting length 19mm, blade length 28mm, overall length 112mm



0103268 **SS**

Werbs DCR Scissors, 5mm angled shaft with rounded tips. Max cutting length 23mm, blade length 30mm, overall length 115mm



Utility

0103227 **SS**

Scissors, straight with sharp tips. Max cutting length 23mm, blade length 30mm, overall length 115mm



0103224 **SS**

Scissors, curved with sharp tips. Max cutting length 23mm, blade length 30mm, overall length 115mm



Bow / Ribbon Style Scissors



0103228 SS

Scissors, straight, short shank with sharp tips.
Max cutting length 18mm, blade length 27mm,
overall length 90mm



0103229 SS

Scissors, curved, short shank with sharp tips.
Max cutting length 18mm, blade length 27mm,
overall length 90mm



0103205 SS

Scissors, straight with sharp tips.
Max cutting length 20mm, blade length 29mm,
overall length 110mm



0103206 SS

Scissors, curved with sharp tips.
Max cutting length 20mm, blade length 29mm,
overall length 110mm



0103202 SS

Scissors, curved with rounded tips. Max
cutting length 24mm, blade length 31mm,
overall length 115mm



0103204 SS

Scissors, straight extra long with rounded tips.
Max cutting length 24mm, blade length 31mm,
overall length 127mm



0103203 SS

Scissors, curved, extra long & strong with rounded
tips. Max cutting length 24mm, blade length
31mm, overall length 127mm

Bow / Ribbon Style Scissors

Utility

0103270 **SS**

Weiss Utility Scissors, straight with sharp tips.
Max cutting length 16mm, blade length 18mm,
overall length 95mm



0103269 **SS**

Weiss Utility Scissors, curved with sharp tips.
Max cutting length 16mm, blade length 18mm,
overall length 95mm



0103271 **SS**

Scissors, Straight, large ring with sharp tips.
Max cutting length 19mm, blade length 28mm,
overall length 100mm



0103272 **SS**

Scissors, Curved, large ring with sharp tips.
Max cutting length 19mm, blade length 28mm,
overall length 98mm



0103275 **SS**

Scissors, Ribbon style, straight with rounded tips.
Max cutting length 18mm, blade length 26mm,
overall length 103mm



0103276 **SS**

Scissors, Ribbon style, curved with rounded tips.
Max cutting length 18mm, blade length 26mm,
overall length 102mm

