

# Beaver-Visitec International

Ophthalmic and Microsurgical Product Catalog



## Wet-Field® Bipolar Instruments

### Wet-Field® Eraser®

Straight Eraser® styles are available in 18G, 20G, 20-23G tapered and 25G designs. Additional specialty products include hollow-core, passive/active aspiration devices with reflux and extrusion cannula modalities - combined with bipolar diathermy capability. Range also includes angled instruments that are designed to provide better tip visualization, range of motion and target approach versatility.

#### Wet-Field® Hemostatic Eraser® Bipolar 18G

10 per box.

**221250 Wet-Field® Hemostatic Eraser® Bipolar**  
18G, bevel tip, straight.



**221251 Wet-Field® Hemostatic Eraser® Bipolar**  
18G, blunt tip, angled 45 degrees.



#### Wet-Field® Hemostatic Eraser® Bipolar 20G

10 per box.

**221260 Wet-Field® Hemostatic Eraser® Bipolar**  
20G, blunt tip, straight.



**221260135 Wet-Field® Hemostatic Eraser® Bipolar**  
20G, blunt tip, angled 135 degrees.




**221260150 Wet-Field® Hemostatic Eraser® Bipolar**  
20G, blunt tip, angled 150 degrees.



#### Wet-Field® Hemostatic Eraser® Bipolar 20-23G Tapered

5 per box.

**221265 Wet-Field® Hemostatic Eraser® Bipolar**  
20-23G tapered, fine tip, straight.

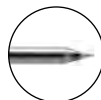



**221265135 Wet-Field® Hemostatic Eraser® Bipolar**  
20-23G tapered, fine tip,  
angled 135 degrees.



**221265150 Wet-Field® Hemostatic Eraser® Bipolar**  
20-23G tapered, fine tip,  
angled 150 degrees.



**221266 Wet-Field® Hemostatic Eraser® Bipolar**  
20-23G tapered, blunt tip, straight.

**22126690 Wet-Field® Hemostatic Eraser® Bipolar**  
20-23G tapered, blunt tip, angled 90 degrees.




**221266135 Wet-Field® Hemostatic Eraser® Bipolar**  
20-23G tapered, blunt tip, angled 135 degrees.






**221266150 Wet-Field® Hemostatic Eraser® Bipolar**  
20-23G tapered, blunt tip, angled 150 degrees.




**Wet-Field® Hemostatic Eraser® Bipolar 23G**

5 per box.



**221268 Wet-Field® Hemostatic Eraser® Bipolar**  
23G, fine tip, straight.

**Wet-Field® Hemostatic Eraser® Bipolar 25G**

Designed to provide fine tip hemostasis without compromising power delivery. 5 per box.

**221267 Wet-Field® Hemostatic Eraser® Bipolar**  
25G, fine tip, straight.

**221267135 Wet-Field® Hemostatic Eraser® Bipolar**  
25G, fine tip, angled 135 degrees.




**221267150 Wet-Field® Hemostatic Eraser® Bipolar**  
25G, fine tip, angled 150 degrees.




## Wet-Field® Bipolar Instruments

### Wet-Field® Hemostatic Eraser® with Aspiration

Designed to reduce the need for instrument exchange during procedure. The Eraser® with active aspirator can be attached to a low vacuum I/A system for variable active aspiration and reflux.

**221261 Wet-Field® Hemostatic Eraser® Bipolar**  
20G, passive aspirator, straight.  
5 per box.



**221262 Wet-Field® Hemostatic Eraser® Bipolar**  
20G, active aspirator with reflux, straight.  
4 per box.



**221262135 Wet-Field® Hemostatic Eraser® Bipolar**  
20G, active aspirator with reflux, angled 135 degrees.  
4 per box.



**221262150 Wet-Field® Hemostatic Eraser® Bipolar**  
20G, active aspirator with reflux, angled 150 degrees.  
4 per box.



### Wet-Field® Coagulators

The Wet-Field® Coagulator is designed to provide precise episcleral, intrascleral, and intraocular hemostasis with reduced peripheral tissue trauma. Controlled Rf diathermy delivers targeted coagulation with the most comprehensive line of bipolar Eraser® instruments, forceps and cables. 1 per box.



Item #	Description
221320	Hemostatic Coagulator (60 Hz), 100-240 volts, with foot switch.
221340	Hemostatic Coagulator, International (50 Hz), 100-240 volts, with foot switch.
232345	Foot Switch. 1 per box.

### Specifications



Weight:	4 lb. 2 oz. (1.9 kg)
Dimensions:	5.5" x 8.5" x 11"
Power Requirements:	100 - 240 VAC, 65 VA, 50/60 Hz
Output Waveform:	445 kilohertz (kHz) square wave variable burst length
Frequency Range:	445 kilohertz
Output Current Range:	36 - 641 milliamperes (mA rms) into a 50 non-inductive load, continuously variable
Output Power:	0.5 - 21 watts into a 50 ohms non-inductive load, continuously variable
Test Load:	50 ohms non-inductive load, neuro-surgical forceps input, using electrosurgical analyzer
Rated Accessory Voltage:	210V Peak



**Wet-Field® Bipolar Forceps**



Bipolar forceps are designed to deliver gentle coagulation with no electrical arcing between the polished instrument tips for minimum coagulum adhesion. The insulated forcep blades are designed to reduce the possibility of electrical shock during operation.



Six available styles provide application versatility. The reliable tip position and consistent operating tension provide quality performance with each use. Reusable. 1 per box.



- 221204 Standard Forceps**  
 Fine point, 15 degree angled tips, handle length 3.5 inches. Designed for pinpoint coagulation.





- 221206 Jewelers Forceps**  
 Fine point, straight tips, handle length 3.5 inches. MacPherson-style fine point microsurgical forceps designed for pinpoint coagulation.



- 221207 Jewelers Forceps**  
 Fine point, 4 mm curved tips, angled 135 degrees, handle length 3.5 inches. MacPherson-style fine point microsurgical forceps.




- 221208 Adson Forceps**  
 Broad tips, handle length 4.2 inches. Designed for facial plastic surgery.



- 221209 Coaptation Forceps**  
 3.5 mm curved tips, angled 118 degrees, handle length 3.5 inches. Designed for sutureless closure of the conjunctiva.

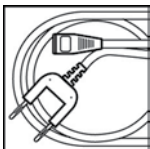

- 221211 Tenzel Forceps**  
 Fine point, 15 degree angled tips, with adjustable tip spacer and built-in stop to prevent tips from touching. Handle length 3.5 inches.

**Wet-Field® Bipolar Cables**



Item #	Description
221202	Wet-Field Cables - Reusable. 3 per box.
22135010	Wet-Field Cables - Sterile, single-use. 10 per box.
22135050	Wet-Field Cables - Sterile, single-use. 50 per box.

Item #	Description
221220	Wet-Field Cables - Reusable, fixed lead. 1 inch distance. 1 per box.

## Accu-Temp® Cautery Devices

### Accu-Temp® Cautery Devices

Accu-Temp® Cauteries provide a broad selection of tip styles for cauterization in a wide variety of ophthalmic procedures.

#### Accu-Temp® High Temperature Cautery

Sterile, single-use. 10 per box.

**8442000 Accu-Temp® High Temperature Cautery**  
1205°C (2200°F),  
1/2 in. shaft with fine tip.



**8443000 Accu-Temp® High Temperature Cautery**  
1205°C (2200°F),  
1/2 in. shaft with loop tip.



**8444000 Accu-Temp® High Temperature Cautery**  
1205°C (2200°F),  
2 in. shaft with loop tip.



**8445000 Accu-Temp® High Temperature Cautery**  
1205°C (2200°F),  
2 in. shaft with fine tip.



**8441500 Accu-Temp® High Temperature Cautery**  
982°C (1800°F),  
1/2 in. shaft with elongated Vasector fine tip.



#### Accu-Temp® Low Temperature Cautery

Sterile, single-use. 10 per box.

**8441000 Accu-Temp® Low Temperature Cautery**  
677°C (1250°F),  
1/2 in. shaft with fine tip.



At Beaver-Visitec International, our promise is to provide health care professionals and service providers worldwide with solutions that exceed expectations.

We are dedicated to providing the best possible customer service and support to our customers. Our team of service, sales and technical support professionals will work with you and your staff, building strong partnerships to ensure your requirements are fully met. We are committed to providing you with the necessary resources, tools, products, and services to enhance your procedures and minimize facility downtime.

Ordering Beaver-Visitec product is easy. Contact a Customer Service Representative via telephone, fax, or email.

### North and South America

Contact our Customer Service and Technical Support teams anytime during normal business hours:  
Monday through Friday, 8:00 am - 6:00 pm Eastern Standard Time

#### Beaver-Visitec International, Inc.

411 Waverley Oaks Road  
Waltham, MA 02452  
USA

Toll-free: 1.866.906.8080  
1.781.906.8080  
Fax: 1.866.906.4304  
Email: [customersupport@beaver-visitec.com](mailto:customersupport@beaver-visitec.com)

### European Countries

Contact our Customer Service and Technical Support teams anytime during normal business hours:  
Monday through Friday, 8:00 am - 5:00 pm Greenwich Mean Time

#### Beaver-Visitec International, Ltd

Centurion Court  
85c Milton Park  
Abingdon, Oxfordshire  
OX14 4RY  
UK

#### Italy

Tel: 39.02.89877520 Option 3  
Fax: 39.02.44386868  
Email: [ordini@beaver-visitec.com](mailto:ordini@beaver-visitec.com)

#### Netherlands

Tel: 31.20808.3202 Option 3  
Fax: 31.20203.1293  
Email: [bestellingen@beaver-visitec.com](mailto:bestellingen@beaver-visitec.com)

#### Belgium

Tel: 32.52.280518 Option 3  
Fax: 32.52.890243  
Email: [bestellingen@beaver-visitec.com](mailto:bestellingen@beaver-visitec.com)

#### Spain

Tel: 34.91.1230335 Option 3  
Fax: 34.91.1849957  
Email: [ordenes@beaver-visitec.com](mailto:ordenes@beaver-visitec.com)

#### France, Luxembourg

Tel: 33.4.85880033 Option 3  
Fax: 33.4.88714057  
Email: [commandes@beaver-visitec.com](mailto:commandes@beaver-visitec.com)

#### United Kingdom, Denmark, Nordic Region

Tel: 44.1865.601.256 Option 3  
Fax: 44.1865.595.761  
Email: [ordersuk@beaver-visitec.com](mailto:ordersuk@beaver-visitec.com)

#### Germany, Switzerland, Austria

Tel: 49.622.16479854 Option 3  
Fax: 49.622.13929883  
Email: [bestellungen@beaver-visitec.com](mailto:bestellungen@beaver-visitec.com)

### Japan and All Other Regions

#### Japan

Tel: 03.5651.2860  
Fax: 03.3662.2242  
Email: [az-center@beaver-visitec.co.jp](mailto:az-center@beaver-visitec.co.jp)

#### All Other Regions

Tel: 44.1865.601.256 Option 3  
Fax: 44.1865.595.761  
Email: [ordersrow@beaver-visitec.com](mailto:ordersrow@beaver-visitec.com)

Beaver-Visitec® product hubs comply with ISO 6009 Color coding and ISO 594 Conical fittings with a 6% luer taper to assure use with syringes or devices.

COLOR	HUB	GAUGE	ID (IN)	ID (MM)	OD (IN)	OD (MM)
Natural* (Clear)		30G	.006	.15	.012	.30
Grey		27G	.008	.20	.016	.40
Brown		26G	.010	.25	.018	.45
Orange		25G	.010	.25	.020	.50
Purple		24G	.012	.30	.022	.55
Blue		23G	.013	.33	.025	.60
Black		22G	.016	.41	.028	.70
Green		21G	.020	.51	.032	.80
Yellow		20G	.024	.60	.036	.90
Cream		19G	.028	.70	.043	1.10

\*Not an ISO color code.





**Beaver-Visitec International, Inc.**

411 Waverley Oaks Road  
Waltham, MA 02452 USA  
Tel: 866.906.8080  
Fax: 866.906.4304

**Beaver-Visitec International Japan, K.K.**

2-13-9 Nihonbashi-Ningyocho  
Chuo-ku, Tokyo 103-0013 Japan  
Tel: 81.3.5651.2861  
Fax: 81.3.3662.2565

**Beaver-Visitec International, Ltd.**

Centurion Court  
85c Milton Park  
Abingdon, Oxfordshire  
OX14 4RY, UK  
Tel: 44.1865.601.265 Option 3  
Fax: 44.1865.595.761